



Quick Reference Card

Cytek® Guava® easyCyte™ Maintenance Guide

For the Guava® easyCyte™ and easyCyte™ HT systems

For Research Use Only. Not for use in diagnostic procedures.

Use this quick start guide when cleaning, running Guava easyCheck, and replacing the flow cell for your Guava easyCyte or easyCyte HT System. For detailed instructions on using the system, refer to the appropriate Guava easyCyte System User's Guide.

Running Guava® Clean

Run Guava® Clean at the start and end of the day, between assays if needed, and to prime the fluid system if you suspect air in the fluid lines. If running lysed whole blood, run two Guava Clean cycles.

Guava® easyCyte™ Systems

1. Click **Cleaning** from the main menu.
2. Click **Start Cleaning**.
3. Place a tube containing 10% bleach in Guava® ICF and a tube of DI water in the sample holder.
4. Load the tube of DI water and click **OK**.
5. When prompted, load the tube of bleach in ICF and click **OK**.
6. When prompted, load a tube of fresh water and click **OK**.

Guava® easyCyte™ HT Systems

1. Click **Cleaning** from the main menu.
2. Click **Start Clean**.
3. Load the following tubes:
 - Seven tubes filled with DI water in positions w1, w2, w3, w4, w5, w6, and tube position 1
 - 100 µL of 100% bleach in tube position 2
 - Two tubes of Guava® ICF in tube positions 3 and 4
 - Tube of DI water in any tube position 5–9 for capillary shutdown
4. Click to select the capillary shutdown position on the plate map, then click **OK**.
5. If you are shutting down the instrument, click **Capillary Shutdown** from the Guava Clean screen. Select a tube from positions 5 through 9 that contains tube of DI water and click **OK**.

Running easyCheck™

Prepare a 1:20 dilution of Guava® easyCheck™ beads. Refer to the Guava easyCheck Kit User's Guide for information. Refer to the appropriate Guava easyCyte™ System User's Guide for information on easyCheck results.

Guava® easyCyte™ Systems

1. Click **easyCheck** from the main menu.
2. Enter the Bead Lot #, Bead Expiration date, and Expected Particles/mL.
3. Mix the easyCheck™ bead sample and load it on the instrument. Click **Run 1st Replicate**.
4. Remove the sample, mix it, then reload it. Click **Run 2nd Replicate**.
5. Remove the sample, mix it again, then reload it. Click **Run 3rd Replicate**.


Guava® easyCyte™ HT Systems

1. Click **easyCheck** from the main menu.
2. Enter the Bead Lot #, Bead Expiration date, and Expected Particles/mL.
3. Click **Start**.
4. Load the following tubes:
 - DI water in positions w2, w4, and w5
 - Empty tubes in positions w3 and w6
 - 100 µL of 100% bleach in position w1
 - easyCheck™ bead sample in any tube or well position
5. Click to select the tube location containing the beads, then click **OK**.


Replacing the Flow Cell

When replacing the flow cell, handle it with care. The capillary is fragile; avoid touching it unnecessarily. Do not force the flow cell into the receptacle.

1. Remove the metal plate from the top of the instrument.
2. Remove the minstac tubing from the clamp and disconnect it from the instrument.

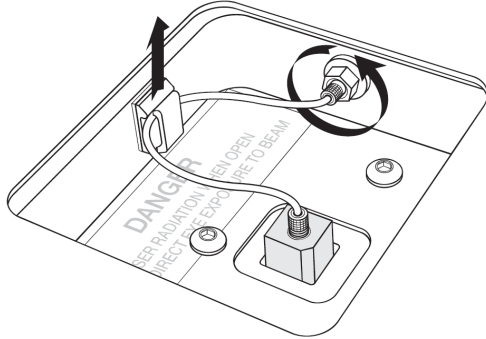


3. Remove the flow cell by grasping it at the top with your fingers and pulling straight up. Do not pull up on the tubing. You can also use the flow cell removal tool to remove the flow cell.
 - Grasp top of the flow cell with the removal tool by placing the clamps into the cutouts at the top of the flow cell, as shown.
 - Pull the tool straight up to remove the flow cell.

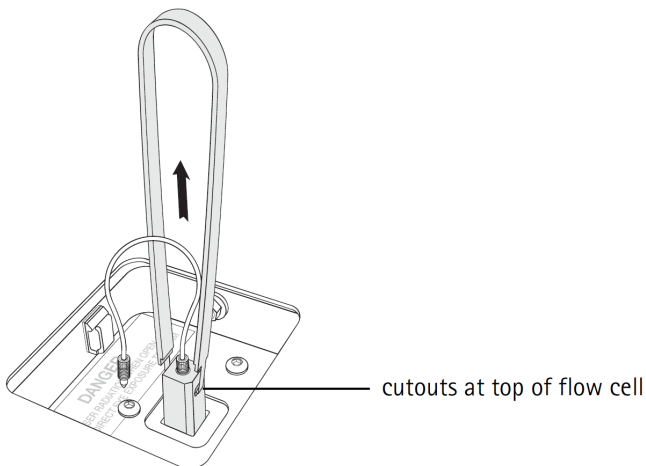


4. Install the new flow cell by correctly positioning it above the instrument and carefully lowering it into the receptacle. It fits one way only—note the beveled edge. Avoid bumping the glass capillary against the instrument or sides of the receptacle. Push down firmly on the top of the flow cell on either side of the

tubing using your fingers until the flow cell clicks into place. Do not press down on the tubing.

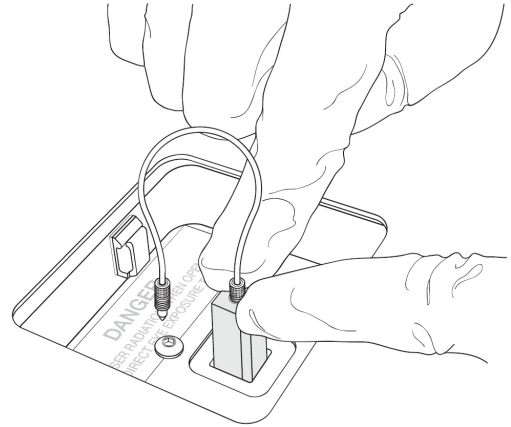

5. Connect the tubing to the instrument, then insert the tubing into the clamp.

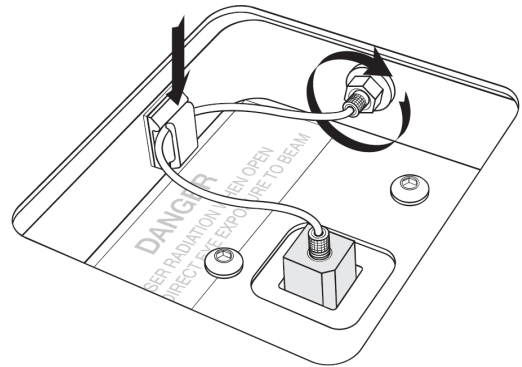

6. Run Guava® Clean to prime the instrument, followed by easy-Check™.

Technical Support

Contact Cytek Biosciences Technical Support by telephone by calling:

U.S. and Other Countries: 1 (877) 922-9835

Europe: +31 (0) 20 765 34 40

Email: technicalsupport@cytekbio.com