



STARTUP

- 1) Turn on Instrument (power switch on back right above cords)
- 2) Open Guavasoft
- 3) Go to Main Menu
- 4) Under Essential Tools, select Cleaning





- 5) Click Start Clean (the tray will eject)
- 6) Fill tubes as described to right
- 7) Select 9 as the Capillary Shutdown position

PS: Ensure that the waste vial is empty and the wash vial is at least half filled with ICF

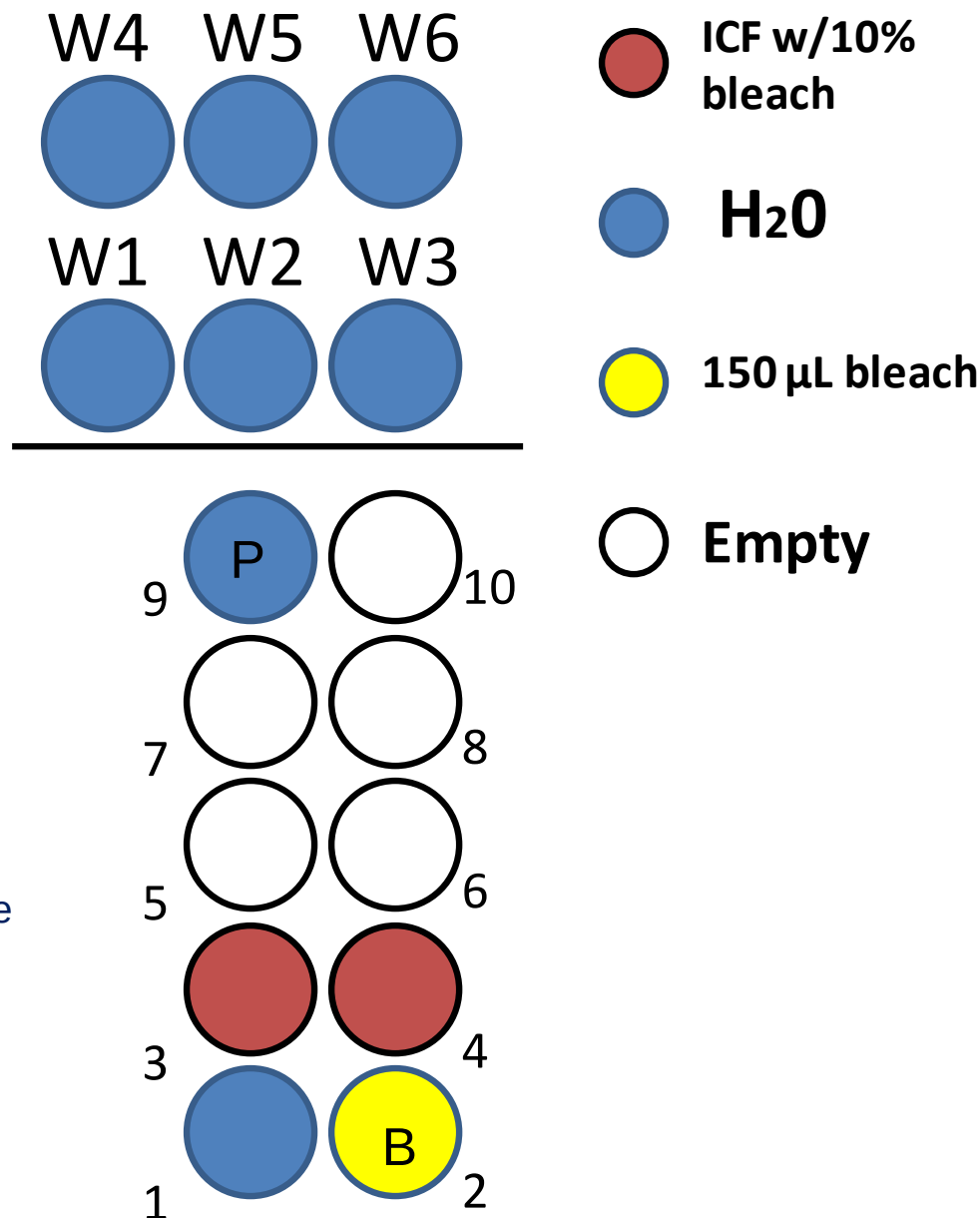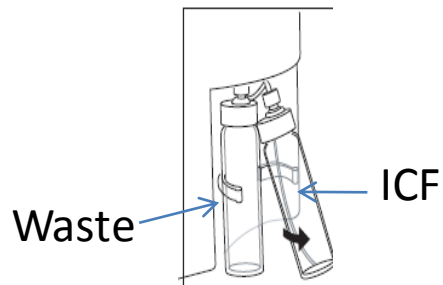

EASYCHECK AND ACQUISITION

- 1) Open Guavasoft
- 2) Open easyCheck
OR
- 3) Open InCyte
- 4) Fill Wash tubes as described on right
- 5) Remember to indicate Capillary Park position on InCyte worklist

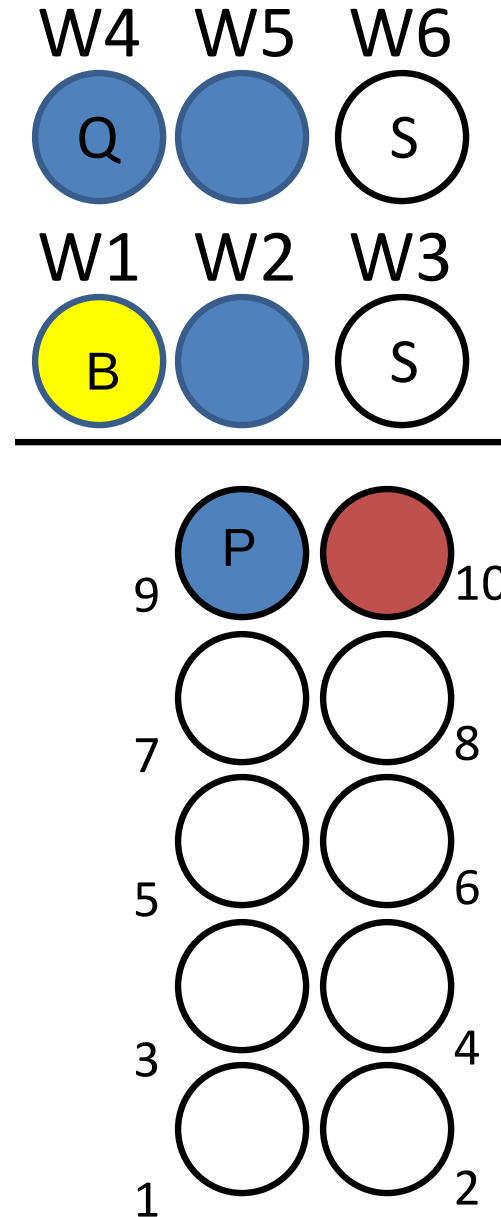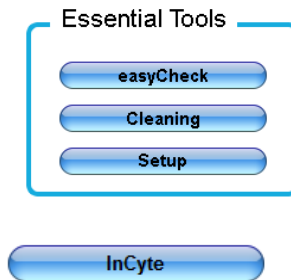

- ICF w/10% bleach
- H2O
- 150 uL bleach
- Empty Tubes for Spin Dry
- Empty (tubes not required)

SHUTDOWN

- 1) Go to Main Menu
- 2) Under Essential Tools, select Cleaning 
- 3) Click Start Clean (the tray will eject)
- 4) Prepare **new** CLEAN tubes as described to right
- 5) Select 9 as the Capillary Shutdown position

